

MedRemind

Your Smart Medicine Reminder

User Guide



Contents

1	Introduction	3
2	Getting Started	3
2.1	Subscribe and Sign In	3
2.2	Create Your Account	5
2.3	Permission Onboarding	6
3	Home Screen	8
4	Medicine Cabinet	11
4.1	Adding a Medicine	11
4.2	Creating a Dose Group	13
5	The Alarm	16
6	Daily Planner	18
7	History	19
8	Settings	20
8.1	Preferences	21
8.2	My Profile	21
9	Tips for Reliable Reminders	22

1 Introduction

MedRemind is a smart medicine reminder app built for patients and caregivers in Bangladesh. It helps you:

- Keep a digital **Medicine Cabinet** of everything you take.
- Organize doses into **Dose Groups** (Morning, Afternoon, Evening, Night) with exact alarm times.
- Get a **full-screen alarm** that rings even when your phone screen is off or you're using another app.
- Track your **adherence** over time with a daily planner and a colour-coded history heatmap.

The app is subscribed to via your mobile operator's billing (Robi/Airtel, through BDApps) at a small daily rate, and authenticates you with your phone number and a password you choose.

2 Getting Started

2.1 Subscribe and Sign In

When you first open MedRemind, you're shown the subscription screen, which explains the price and the three headline features.

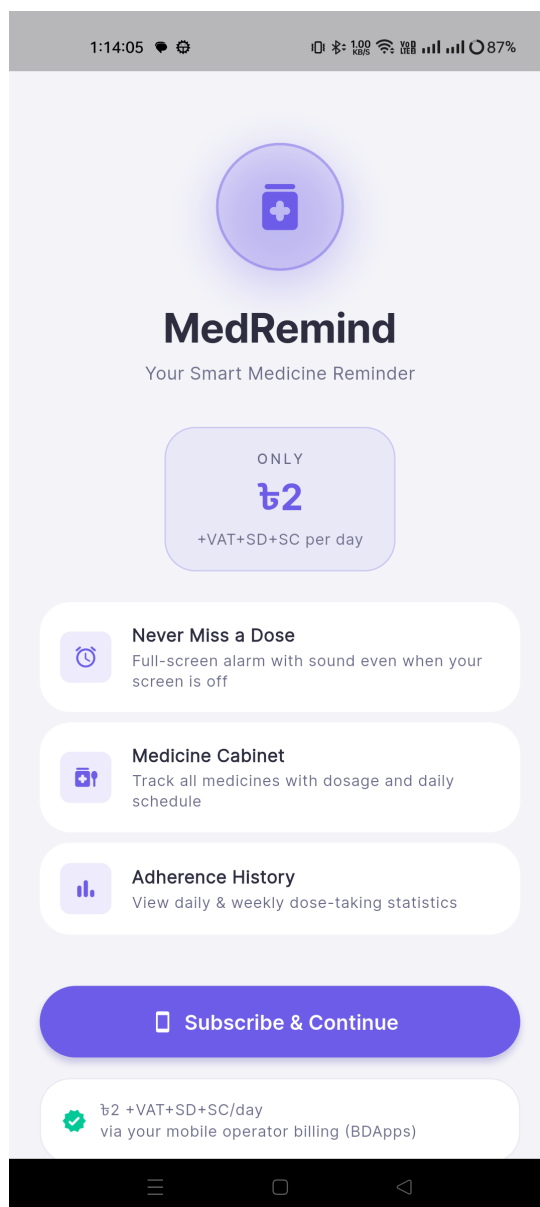


Figure 1: Subscription / welcome screen

Tapping **Subscribe & Continue** takes you to the sign-in screen, where you enter your mobile number. If the number is new, you'll be taken to registration; if it already has an account, you'll be asked for your password instead.


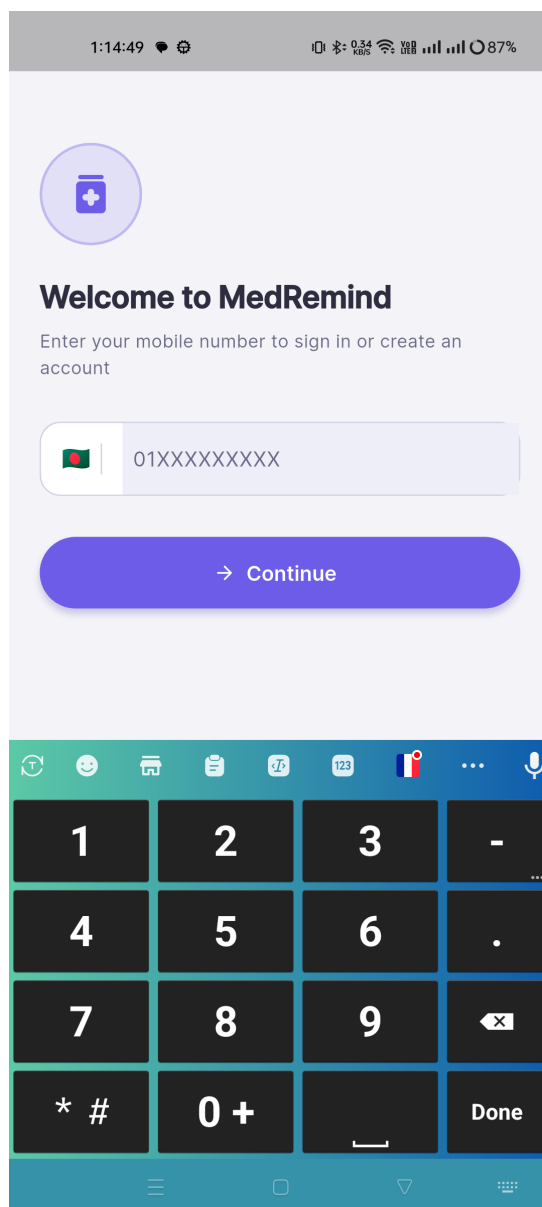


Figure 2: Enter your mobile number to sign in or register

2.2 Create Your Account

New numbers are asked for a name and a password (minimum 6 characters).

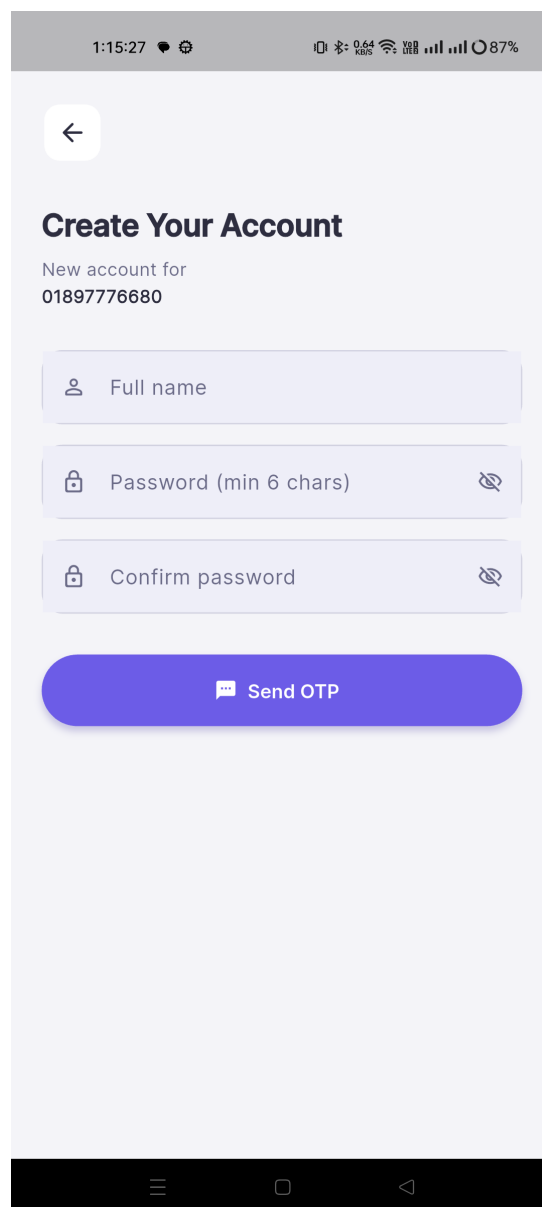

1:15:27 0.64 KB/s 87%

←

Create Your Account

New account for
01897776680

Full name

Password (min 6 chars)

Confirm password

Send OTP

Figure 3: Create Your Account

After submitting, MedRemind sends a one-time SMS code to verify you own the number before your account is created. Enter the 6-digit code on the OTP screen to finish registration.

2.3 Permission Onboarding

A short three-screen walkthrough introduces the app's core value before you grant any permissions:

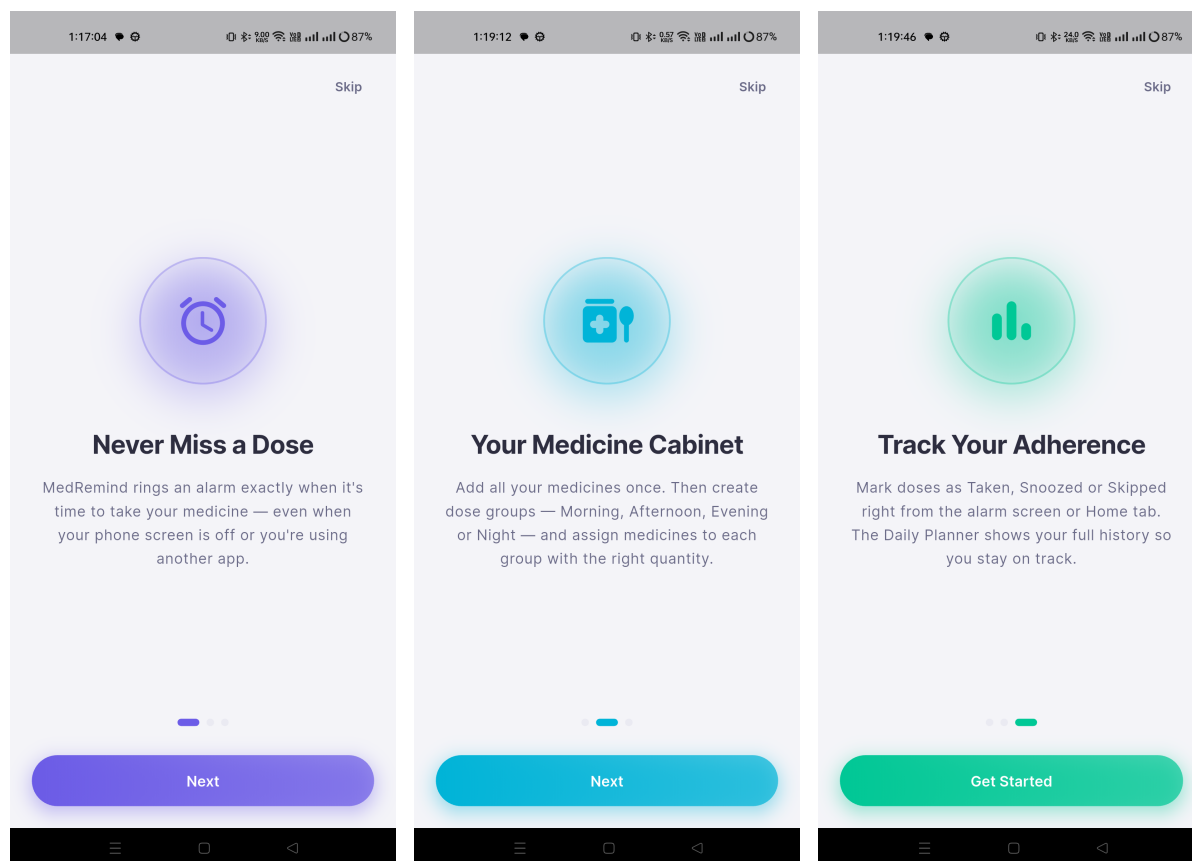


Figure 4: Onboarding: Never Miss a Dose, Your Medicine Cabinet, Track Your Adherence

Right after, MedRemind asks for the permissions it needs to actually ring alarms reliably: notifications, exact alarms, and an exemption from battery optimization (so the OS doesn't kill the alarm in the background).

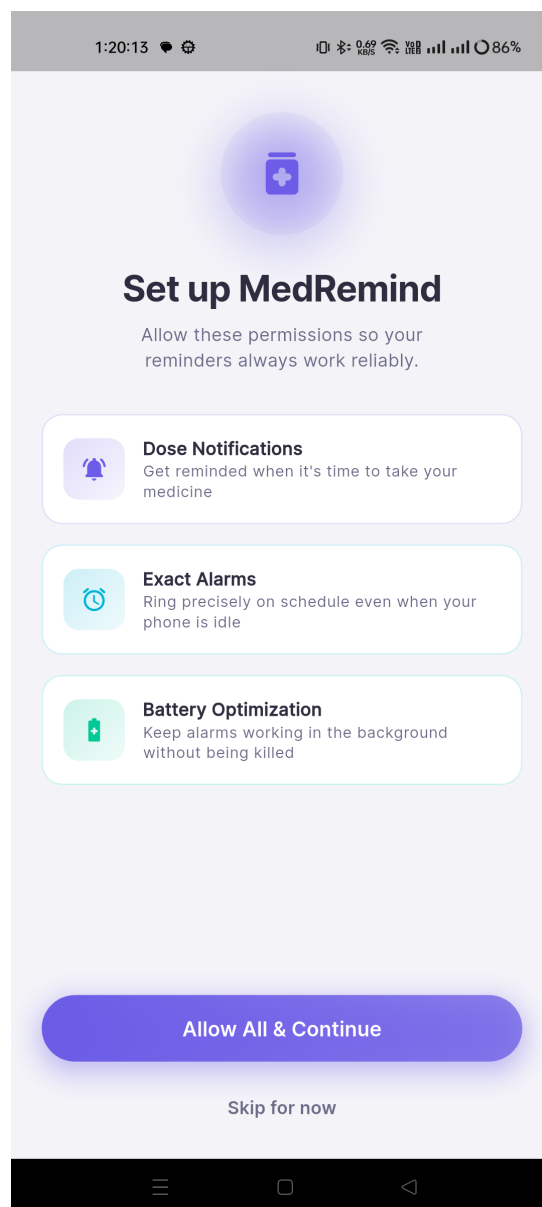


Figure 5: Set up MedRemind — permission requests

You can tap **Allow All & Continue** to grant everything at once, or **Skip for now** and grant permissions later from Settings. Skipping these, especially *Exact Alarms* and *Battery Optimization*, may cause reminders to ring late or not at all.

3 Home Screen

The Home tab is what you land on after signing in. It shows:

- A **Today's Pills** summary card with a progress ring.
- Filter chips: *All / Morning / Afternoon / Evening / Night*.
- A 7-day date strip to jump to a different day.
- The list of doses due that day, each with **Taken**, **Snooze**, and **Skip** actions.

Before any medicines are added, the Home screen shows an empty state prompting you to add your first medicine.

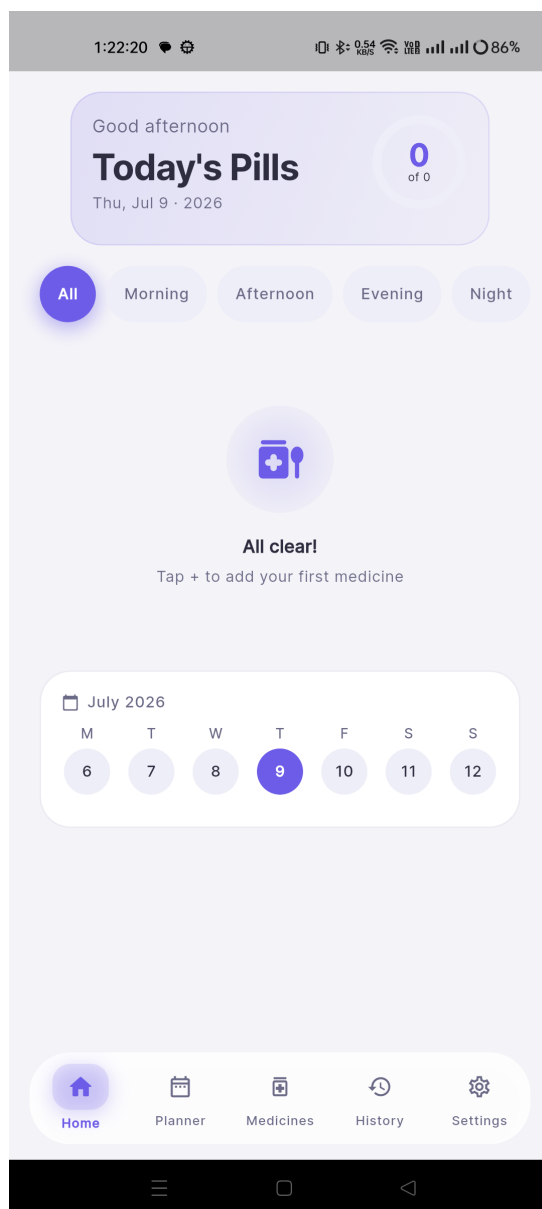


Figure 6: Home tab — empty state

Once a dose group exists, its doses appear as cards on Home. Each dose card lets you mark it Taken, Snooze it for a few minutes, or Skip it outright.

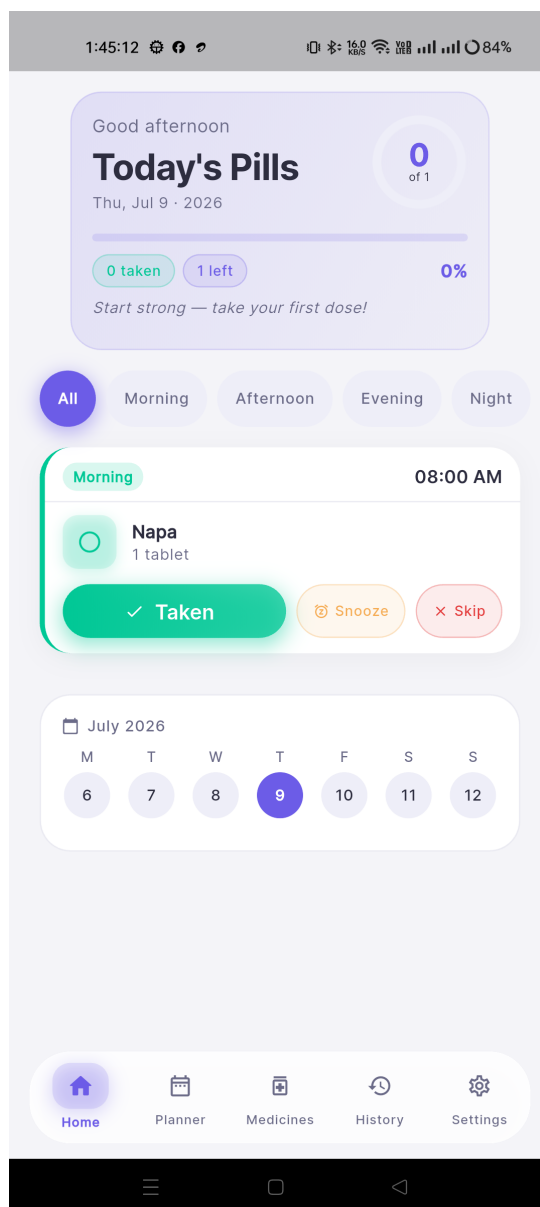


Figure 7: A due dose, with Taken / Snooze / Skip actions

Once every dose for the day is marked Taken, the progress ring fills to 100% and the day's remaining schedule (today's completed doses, and tomorrow's upcoming ones) is shown below.

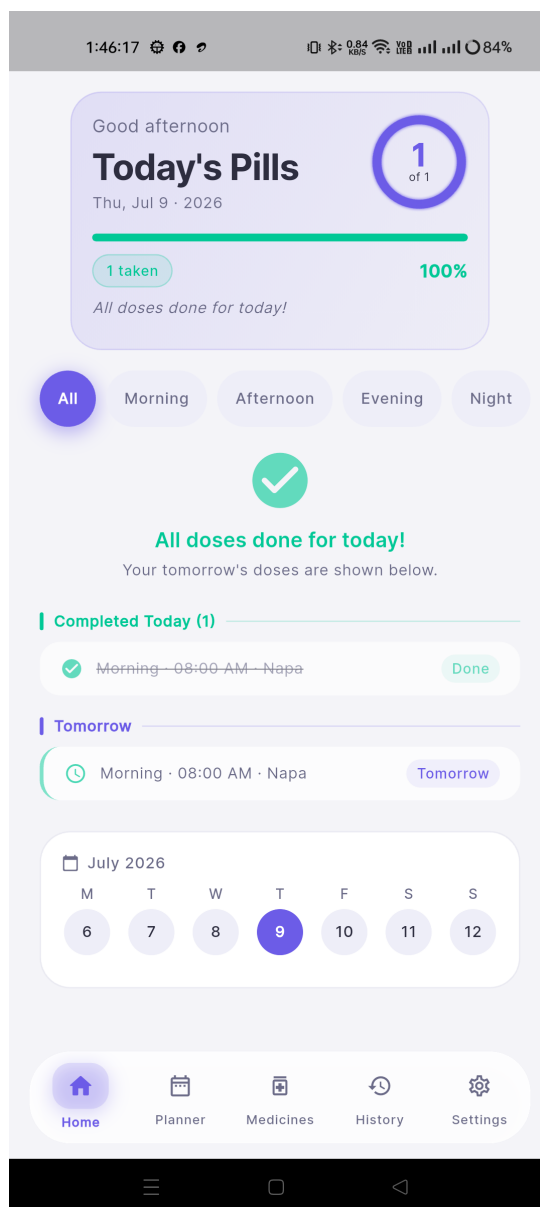


Figure 8: All doses done for today

4 Medicine Cabinet

The **Medicines** tab is where you manage what you take and when. It has two sub-tabs: *Medicines* (the raw list of medicines you take) and *Schedules* (the dose groups that actually trigger alarms).

4.1 Adding a Medicine

Tap + **Add Medicine** from the empty Medicines list to open the add form.

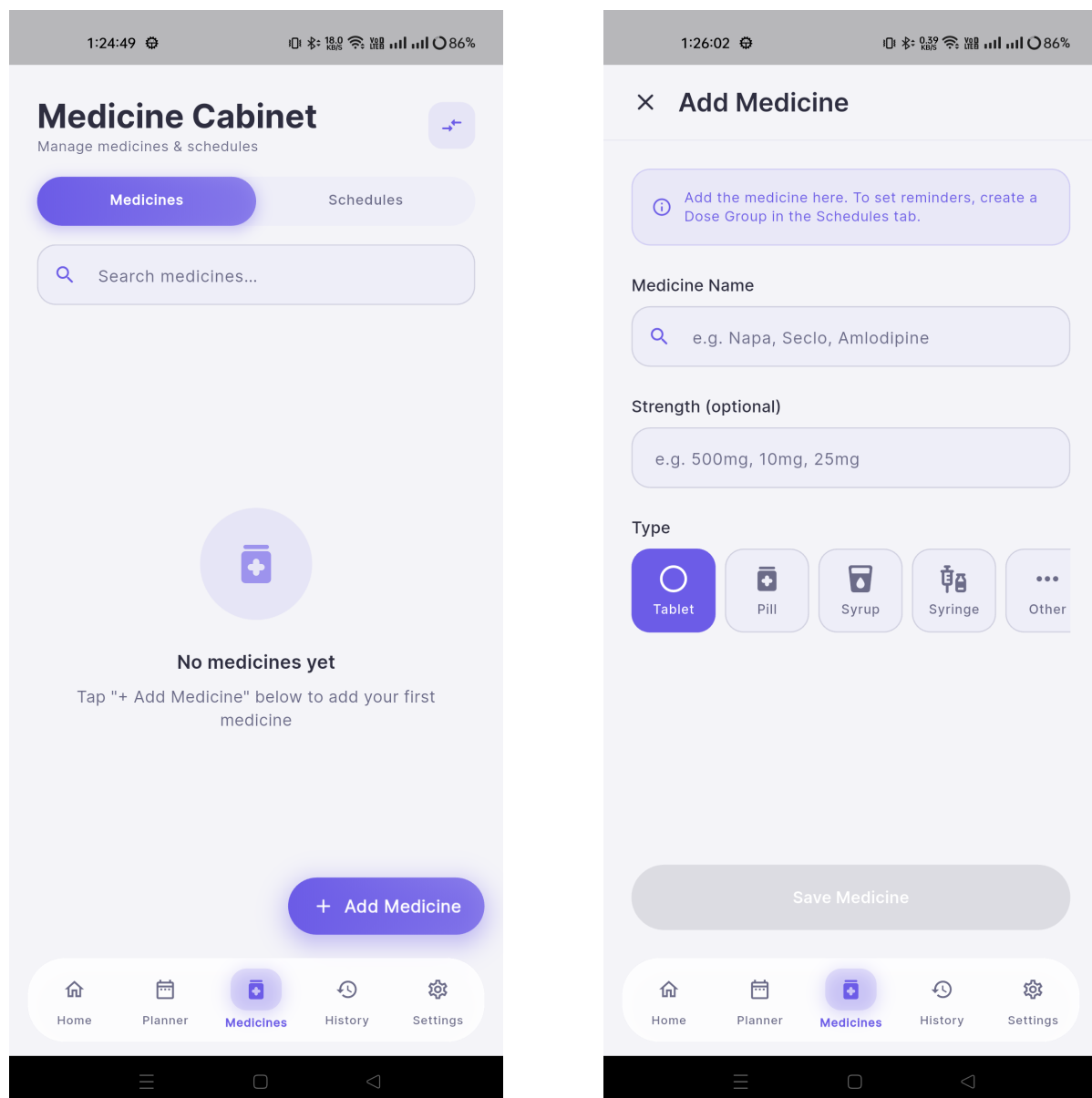


Figure 9: Medicine Cabinet (empty) and the Add Medicine form

As you type a medicine name, MedRemind suggests matches from a built-in database of common Bangladeshi brand names, along with the generic drug name (e.g. *Napa* → Generic: Paracetamol). You can optionally set a strength (e.g. 500mg) and a form (Tablet, Pill, Syrup, Syringe, or Other).

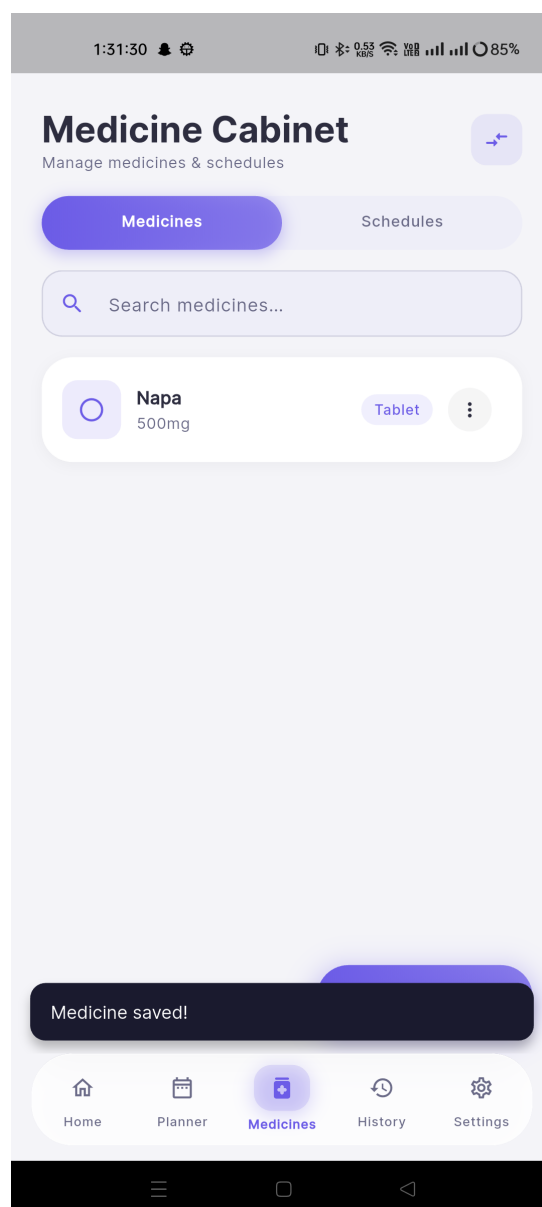


Figure 10: Medicine saved and listed in the cabinet

Note: adding a medicine here does *not* by itself create any reminders. To actually get an alarm, you need to add the medicine to a *Dose Group*, covered next.

4.2 Creating a Dose Group

Switch to the **Schedules** sub-tab and tap **+ New Group**. A dose group bundles together:

- A **Dose Time** label (Morning / Afternoon / Evening / Night) — purely descriptive, shown throughout the app.
- An exact **Alarm Time**.
- An optional **Meal Relation** (No preference / Before meal / After meal) — shown as a note on the reminder.
- **Repeat Days** — every day by default, or specific weekdays.

- One or more **Medicines**, each with its own quantity.

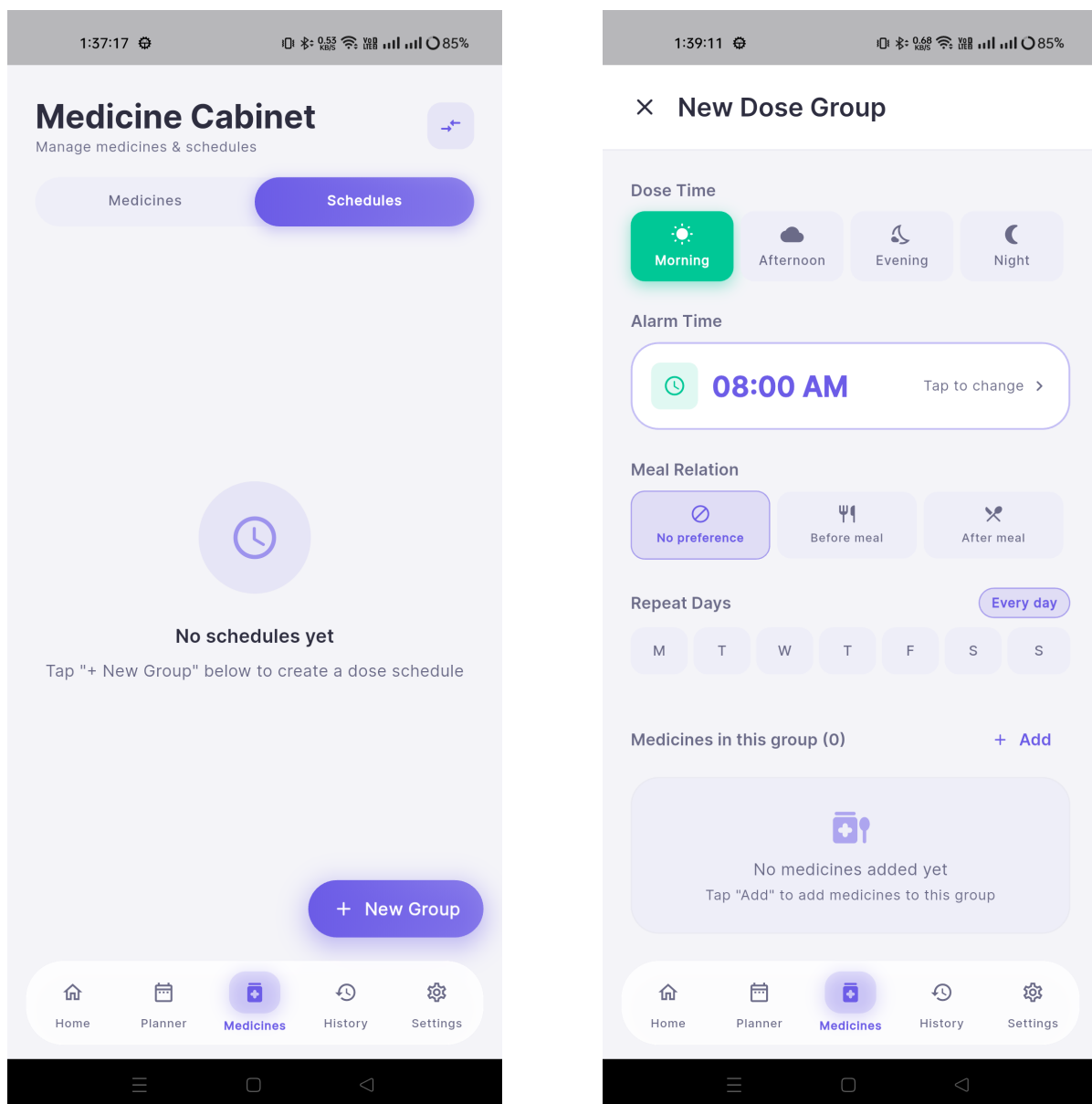


Figure 11: Schedules tab (empty) and the New Dose Group form

Tap **Add** under "Medicines in this group" to pick from your cabinet. Each added medicine gets its own $-/+$ quantity stepper.

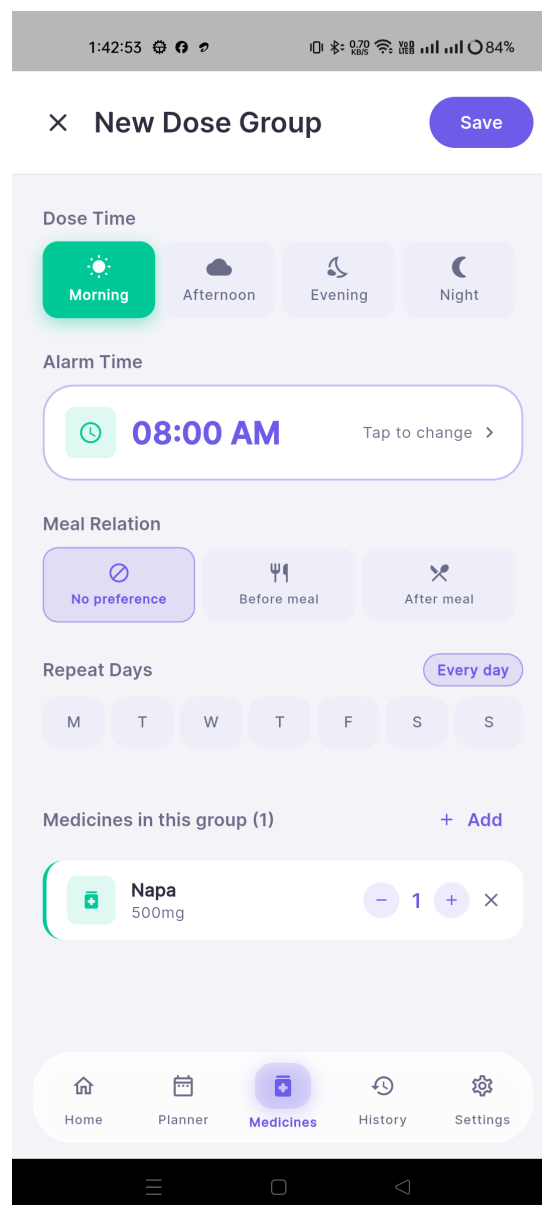


Figure 12: A medicine added to the dose group, with quantity control

Once saved, the group appears in the Schedules list, and its dose will now show up on the Home tab and in the Daily Planner at the time you chose.

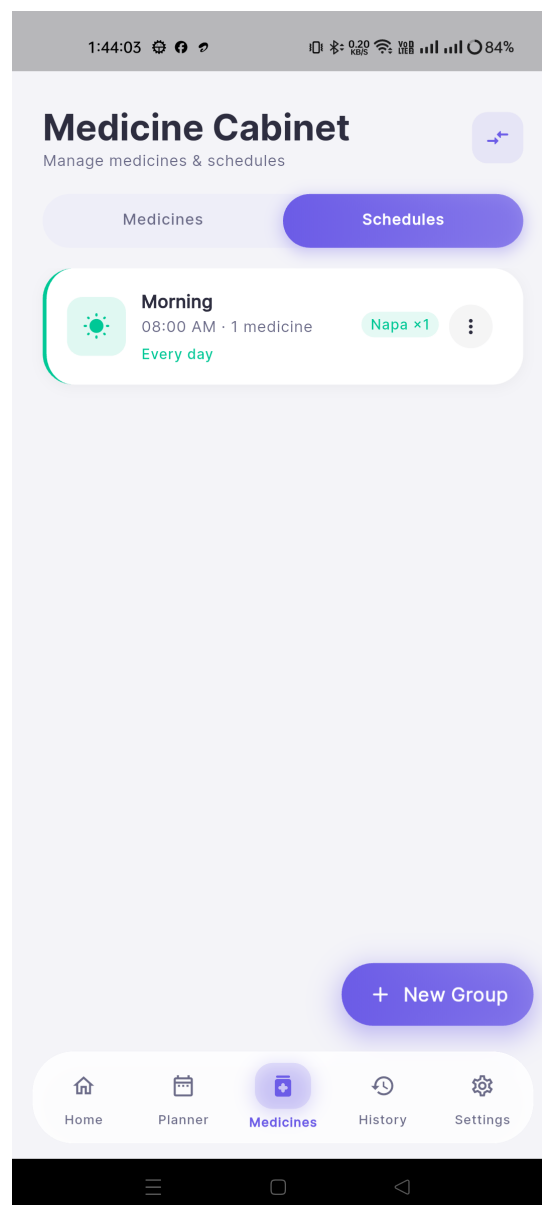


Figure 13: Saved dose group: Morning, 08:00 AM, Napa ×1, every day

5 The Alarm

When a dose group's alarm time arrives, MedRemind shows a full-screen alarm — even over the lock screen or another app — with the medicine name and quantity, and three actions: **Dismiss** — **Taken**, **Snooze 5 min**, and **Skip**. Swiping up also marks it Taken.

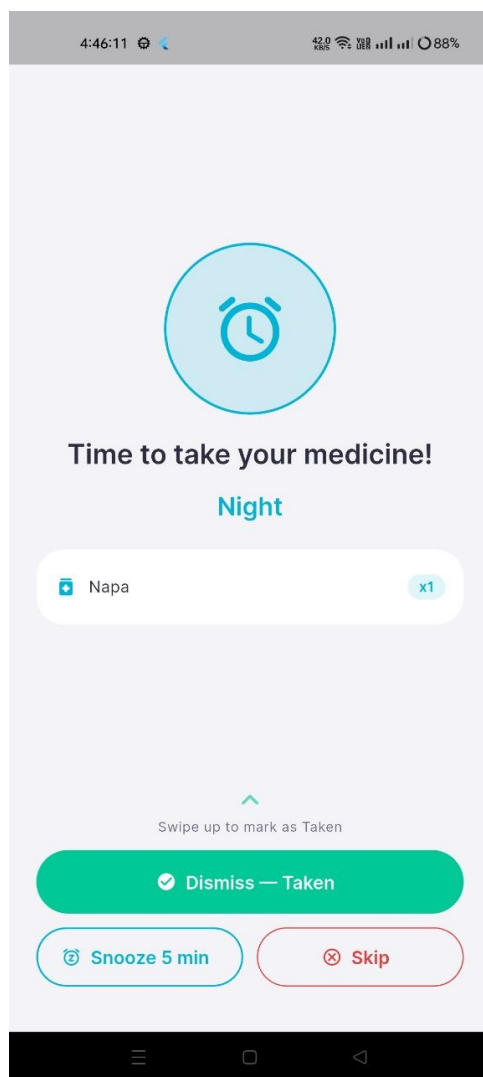


Figure 14: Full-screen medicine alarm

If you don't see or act on the alarm right away, a regular notification is also posted with the same three actions (Dismiss/Taken, Snooze, Skip) available directly from the notification shade.

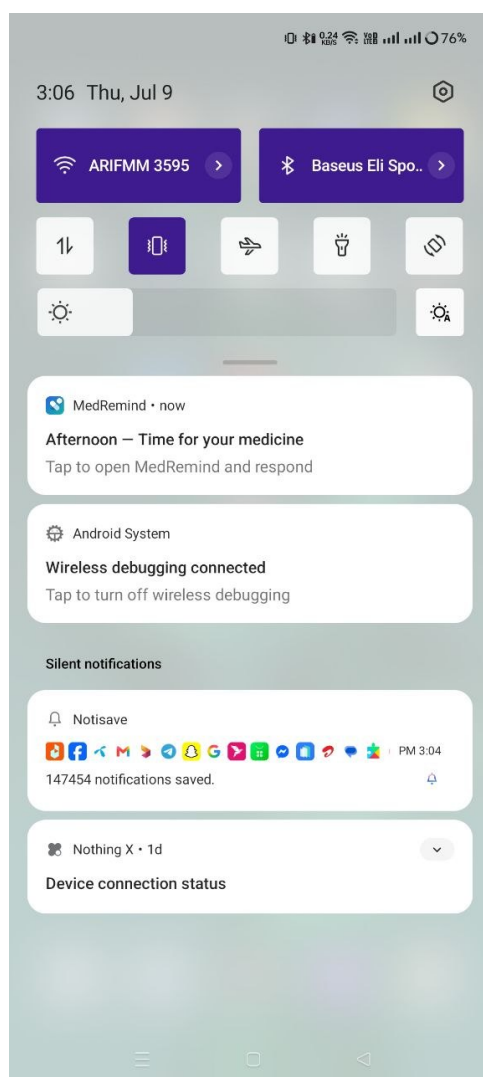


Figure 15: Medicine alarm notification, actionable without opening the app

6 Daily Planner

The **Planner** tab shows a full month calendar. Each day carries a coloured dot summarizing that day's adherence:

- **Green** — All doses taken
- **Orange** — Partial (some doses taken, some not)
- **Red** — Missed
- **Purple** — Scheduled (future day)

Selecting a day shows a summary card and the full list of that day's doses with their status.

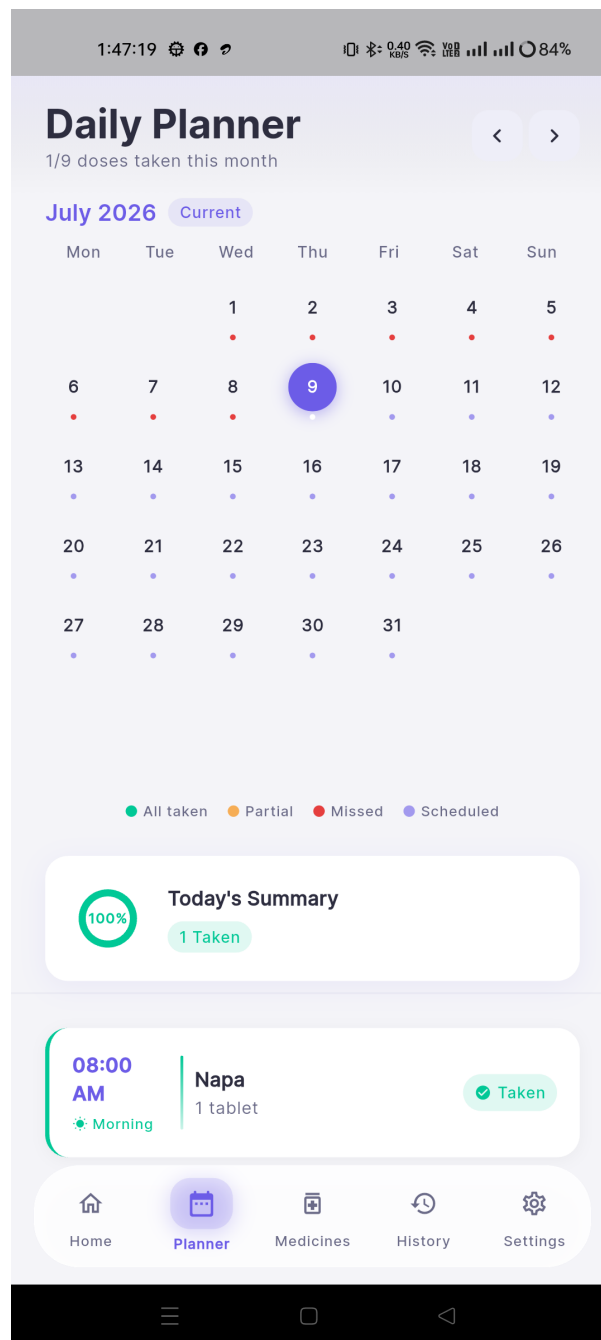


Figure 16: Daily Planner with adherence heatmap

7 History

The **History** tab summarizes your adherence over time: a "this week" percentage with taken/missed/skipped counts, a monthly adherence heatmap (same colour legend as the Planner), and a chronological **Dose log** listing every recorded dose event.

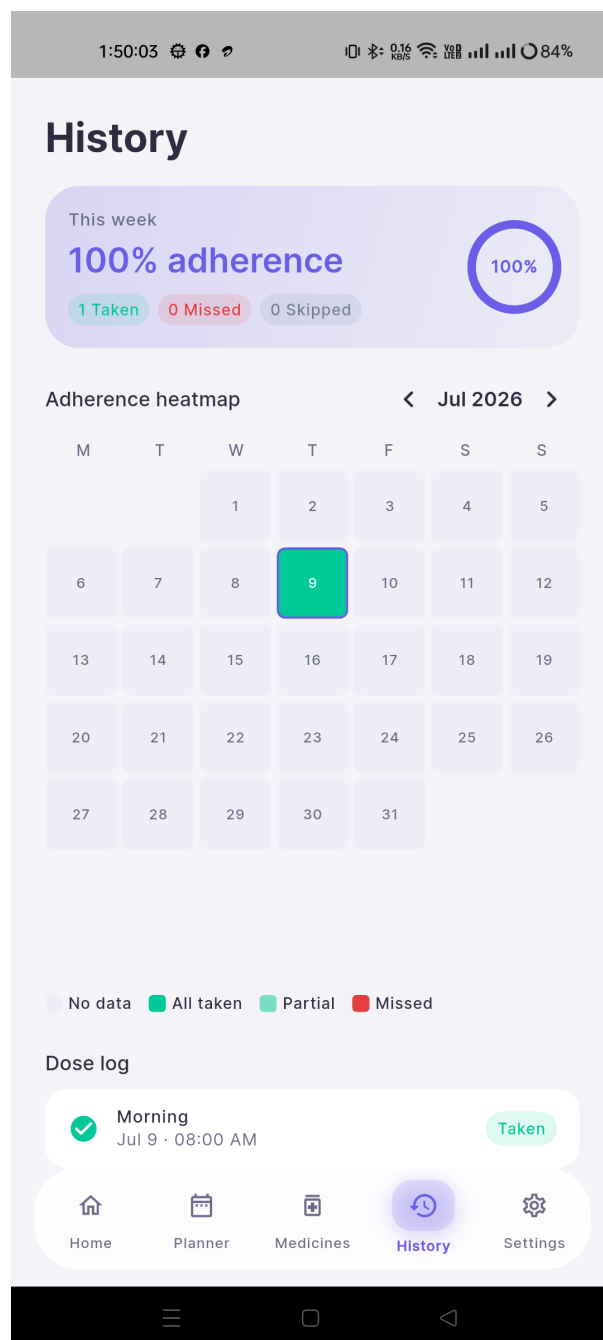


Figure 17: History — adherence heatmap and dose log

8 Settings

The **Settings** tab holds your profile and app preferences.

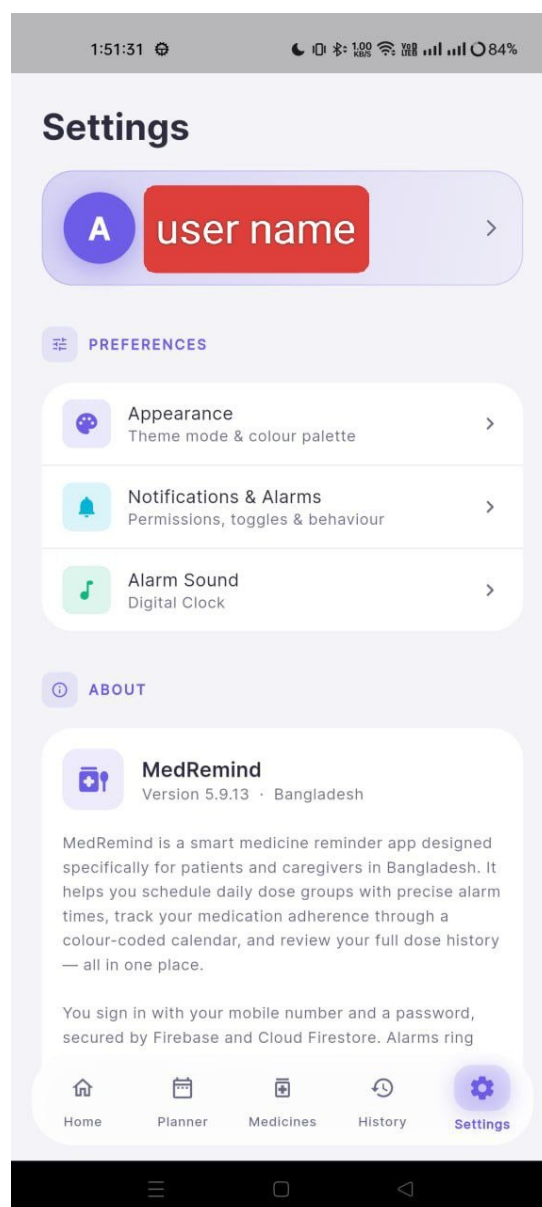


Figure 18: Settings — profile, preferences, and about

8.1 Preferences

- **Appearance** — theme mode and colour palette.
- **Notifications & Alarms** — permissions, toggles, and behaviour for how reminders fire.
- **Alarm Sound** — choose the ringtone used for the full-screen alarm.

8.2 My Profile

Tapping your name at the top of Settings opens your profile, showing your mobile number and subscription status/expiry, plus account actions:

- **Logout** — sign out of this device only; your account and data remain.
- **Change Password** — update your password (requires your current one).

- **Unsubscribe** — cancels your BDApps subscription *and* permanently deletes your account and all its data. This cannot be undone.

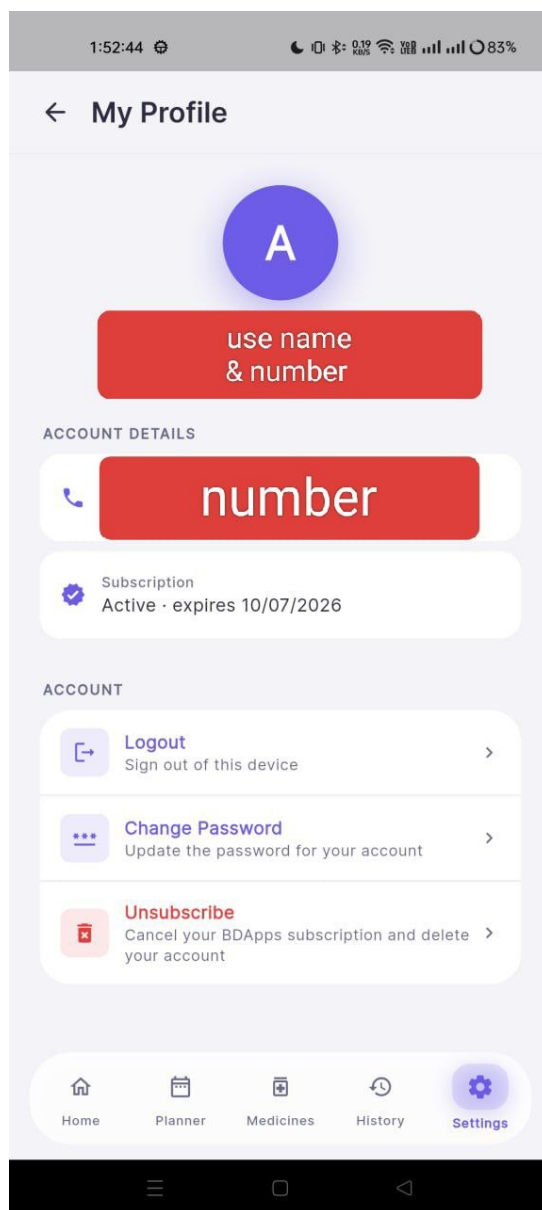


Figure 19: My Profile — account details and account actions

9 Tips for Reliable Reminders

1. Grant **Exact Alarms** and disable **Battery Optimization** for MedRemind when prompted — without these, some Android phones delay or silently drop alarms.
2. Keep the app installed and signed in; alarms are scheduled locally on your device and don't need an active internet connection to ring.
3. Use **Dose Groups** rather than one-off reminders — grouping medicines by time keeps the Home screen and Planner clean and accurate.
4. Check the **History** tab periodically; a string of missed or skipped doses is worth discussing with your doctor or caregiver.